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Instruction for BrainAccess EEG Devices

Installation (Windows 10/11)

1. Download the latest version of [BrainAccess](#) Board for Windows 10/11.
2. Unzip the file and run the BrainAccessBoard_v2.6.1.exe installer located inside.
3. Follow the on-screen instructions.


Installation (Linux)

1. Download the latest version of [BrainAccess](#) Board for Linux.
2. Unzip the file and run the BrainAccessBoard_v2.6.1.deb file located inside.
3. Follow the on-screen instructions.


Connecting the Device

1. Turn on the device by holding the button until the LED lights up blue (the HALO headband will additionally vibrate).
2. Enable Bluetooth on your computer or use the Bluetooth adapter included in the kit (DO NOT PAIR/CONNECT THE DEVICE yourself via the system settings).
3. Launch the installed BrainAccess Board program.

4. Select the device you want to connect to from the list.

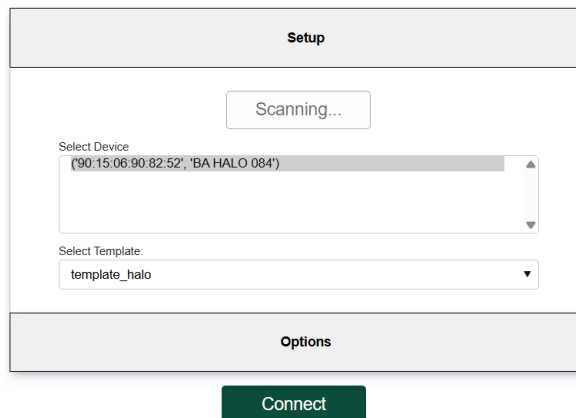5. If the connection is successful, the device should appear on the list.

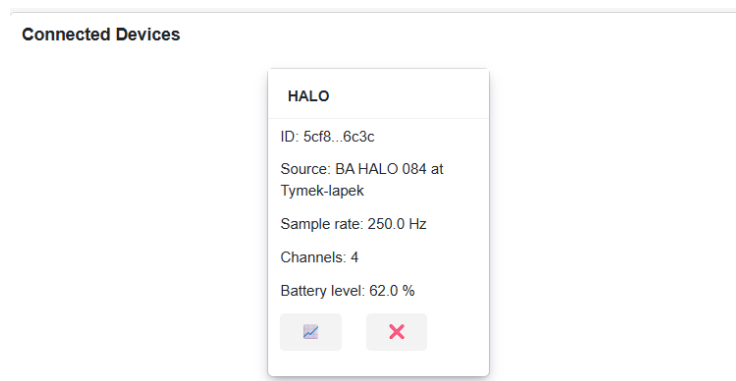Troubleshooting

1. **I can't see my device on the list.**

- Check if your device is turned on and not out of battery.
- Check if Bluetooth is enabled or if you are using the adapter.
- Try closing the BrainAccess Board program, then first turn on the headband and Bluetooth, wait for 10 seconds, and launch the application again.

2. **The signal is very weak.**

- Make sure the electrodes have good contact with the head.
- In the case of the HALO headband, check that no hair is under the flat forehead electrodes.

3. **Feel free to ask our technical stuff for any help with EEG devices!**